

Σ. ΚΟΥΖΟΥΝΑΣ & ΣΙΑ Ε.Ε. - ALCOFARM MEDICAL

STUDY REPORT 2020-7622/20 23 00657

Alco X Quat

**SUSPENSION TEST
ACCORDING TO EN 13624:2013
(Phase 2 step 1)**

Chemical disinfectants and antiseptics
Quantitative suspension test for the evaluation of fungicidal or
yeastocidal activity in the medical area - Test method
and requirements (phase 2, step 1)

SEPTEMBER 2020

STUDY REPORT 2020-7622/20 23 00657

SUSPENSION TEST ACCORDING TO EN 13624:2013

Chemical disinfectants and antiseptics - Quantitative suspension test for the evaluation of yeasticidal activity in the medical area - Test method and requirements
(Phase 2 step 1)

TEST PRODUCT IDENTIFICATION

PRODUCT NAME	:	Alco X Quat
ACTIVE SUBSTANCES	:	Ethyl Alcohol Denat. 59% w/v (73,75% v/v) Benzalkonium Chloride (50%): 0.8%
APPEARANCE OF THE PRODUCT	:	Liquid
STORAGE CONDITIONS	:	Room temperature, darkness
TEST CONDITIONS	:	Test conducted at 20°C ± 1 °C
LOT	:	Not provided
METHOD	:	EN 13624:2013
CONTACT TIME	:	5 minutes
DILUTIONS	:	Undiluted (80%), 50%, 1%
STUDY SPONSOR	:	Σ. KOYZOYNAΣ & ΣΙΑ Ε.Ε. -ALCOFARM MEDICAL
PRODUCT SUPPLIER	:	Σ. KOYZOYNAΣ & ΣΙΑ Ε.Ε. - ALCOFARM MEDICAL
PRODUCT MANUFACTURER	:	Σ. KOYZOYNAΣ & ΣΙΑ Ε.Ε. - ALCOFARM MEDICAL
RECEIPT DATE	:	10/07/2020
STUDY PERIOD	:	11/08/2020-13/08/2020
LAB ID	:	2020-7622/20 23 00657

SCOPE

This European Standard specifies a test method and the minimum requirements for fungicidal or yeasticidal activity of chemical disinfectant and antiseptic products that form a homogeneous, physically stable preparation when diluted with hard water, or - in the case of ready-to-use products - with water. Products can only be tested at a concentration of 80 % or less (97 % with a modified method for special cases) as some dilution is always produced by adding the test organisms and interfering substance.

This European Standard applies to products that are used in the medical area in the fields of hygienic handrub, hygienic handwash, surgical handrub, surgical handwash, instrument disinfection by immersion, and surface disinfection by wiping, spraying, flooding or other means.

This European Standard applies to areas and situations where disinfection or antiseptics is medically indicated. Such indications occur in patient care, for example:

- in hospitals, in community medical facilities and in dental institutions;
- in clinics of schools, of kindergartens and of nursing homes;

and may occur in the workplace and in the home. It may also include services such as laundries and kitchens supplying products directly for the patients.

PRINCIPLE

A sample of the product as delivered and/or diluted with hard water (or water for ready to use products) is added to a test suspension of yeasts in a solution of an interfering substance. The mixture is maintained at the temperature (θ) and for the chosen contact time (t). At the end of this contact time, an aliquot is taken; the yeasticidal action in this portion is immediately neutralized or suppressed by a validated method. The method of choice is dilution-neutralization. If a suitable neutralizer cannot be found, membrane filtration is used. The numbers of surviving yeasts in each sample are determined and the reduction is calculated. Handwash products are always prediluted with hard water. The resulting solution is regarded as a ready-to-use product.

TEST CONDITIONS

1. Product type: Surface Disinfection.
2. The following procedure was performed in water bath at 20 °C.
3. The test product was tested at 5 minutes contact time.
4. Interfering substance: A final concentration of 0.3g/L bovine albumin was used for testing (clean conditions).
5. Neutralization Method used: Dilution neutralization.
6. Neutralizer used: LPT Dilution Broth containing polysorbate 80.
7. According to EN 13624, products shall be tested at a minimum of three different concentrations to include one concentration in the active range and one concentration in the non-active range. In this case the product was tested: Undiluted, 50% 1%.

TEST MICROORGANISMS

Candida albicans

ATCC 10231

Aspergillus brasiliensis

ATCC 16404

FUNGICIDAL ACTIVITY FOR SURFACE DISINFECTION PRODUCTS

The product shall be deemed to have passed the EN 13624 Standard if it demonstrates in a valid test for surface disinfection products at least a 4 lg reduction within max. 5 min (or between 6 min and 60 min for products used on surfaces which do not require an action within 5 min or shorter) at min. 4 °C and max. 30 °C with the chosen interfering substance (clean or dirty conditions) under the conditions defined by this standard when the test organisms are *Aspergillus brasiliensis* and *Candida albicans*.

ASSAY ACCEPTANCE CRITERIA

1. Test Suspension (N) is between 1.5 to 5.0 X 10⁷ CFU per mL (7.17≤log No≤7.70)
2. No (N/10) is between 1.5 to 5.0 X 10⁶ CFU per mL (6.17≤log No≤6.70)
3. Validation Suspension=Nv is between 3.0 x 10² and 1.6 x 10³.
4. Neutralizer control= NvB is between 3.0 x 10⁴ and 1.6 x 10⁵.
5. Nvo (Nv/10) is between 30 and 160.
6. Na is the number of survivors (cells) per ml in the test mixture at the end of contact time.
7. R (log reduction) = No - Na
8. Average recovery values for the experimental conditions control (A) were equal to or greater than 0.5 times the Validation Suspension (Nvo)
9. Average recovery values for the Neutralizer control (B) were equal to or greater than 0.5 times the Validation Suspension (Nvo) or NvB/1000.
10. Average recovery values for the Method Validation control (C) were equal to or greater than 0.5 times the Validation Suspension (Nvo)

TEST RESULTS FOR *Candida albicans* (YEASTICIDAL SUSPENSION TEST)

Test suspension

Test - suspension (N and No)			
N	Vc1	Vc2	x mean 2.87E+07
10 ⁻⁵	298	276	
10 ⁻⁶	27	31	
			log N 7.46
			No (N/10) 2.87E+06
			log No 6.46
			6,17 < = logNo < = 6,70 Yes

Validation and controls

Validation suspension (Nvo)	Experimental conditions (A)	Neutralizer control (B)	Method validation (C) Product conc.:	Undiluted (80%)
VC 1 64 VC 2 57 x mean 60.5	VC 1 61 VC 2 68 x mean 64.5	VC 1 69 VC 2 75 x mean 72	VC 1 74 VC 2 79 x mean 76.5	
30*x mean of Nvo < 160? Yes	x mean of A is > 0,5*x mean of Nvo? Yes	x mean of B is > 0,5*x mean of Nvo or Nvb/1000? Yes	x mean of C is > 0,5*x mean of Nvo? Yes	
Validation suspension (Nvb)				
VC 1 62 VC 2 54 x mean 58				
30*x mean of Nvb < 160? Yes				

Test Results

Product concentration (%)	Contact time	Dilution step	Vc 1	Vc 2	Average of Vc 1 and Vc 2	Na= average x10	log Na	log No	log Reduction (No-Na)	Criteria	Result
Undiluted (80%)	5 min	10 ⁰	0	0	< 14	< 140	< 2.15	6.46	> 4.31	≥4	PASS TEST
		10 ⁻¹	0	0							
50%	5 min	10 ⁰	0	0	< 14	< 140	< 2.15	6.46	> 4.31	≥4	PASS TEST
		10 ⁻¹	0	0							
1%	5 min	10 ⁰	> 330	> 330	< 3300	< 33000	< 4.52	6.46	> 1.94	≥4	FAILS TEST
		10 ⁻¹	> 330	> 330							

TEST RESULTS FOR *Aspergillus brasiliensis* (FUNGICIDAL SUSPENSION TEST)

Test suspension

Test - suspension (N and No)			
N	Vc1	Vc2	x mean
10^{-6}	21	24	2.27E+07
10^{-7}	3	2	log N
			7.36
			No (N/10)
			2.27E+06
			log No
			6.36
			6,17 < = logNo < = 6,70
			Yes

Validation and controls

Validation suspension (Nvo)			Experimental conditions (A)			Neutralizer control (B)			Method validation (C) Product conc.:			Undiluted (80%)
VC 1	49	x mean	VC 1	59	x mean	VC 1	48	x mean	VC 1	44	x mean	
VC 2	54	51.5	VC 2	55	57	VC 2	50	49	VC 2	54	49	
30*x mean of Nvo < 160?			x mean of A is > 0,5*x mean of Nvo?			x mean of B is > 0,5*x mean of Nvo or Nva/1000?			x mean of C is > 0,5*x mean of Nvo?			
Yes			Yes			Yes			Yes			
Validation suspension (Nva)												
VC 1	40	x mean										
VC 2	37	38.5										
30*x mean of Nva < 160?			Yes									

Test Results

Product concentration (%)	Contact time	Dilution step	Vc 1	Vc 2	Average of Vc1 and Vc2	Na= average x10	log Na	log No	log Reduction (No-Na)	Criteria	Result
Undiluted (80%)	5 min	10^0	> 165	> 165	> 1650	> 16500	> 4.22	6.36	< 2.14	≥4	FAILS TEST
		10^{-1}	> 165	> 165	> 1650	> 16500	> 4.22	6.36	< 2.14	≥4	FAILS TEST
50%	5 min	10^0	> 165	> 165	> 1650	> 16500	> 4.22	6.36	< 2.14	≥4	FAILS TEST
		10^{-1}	> 165	> 165	> 1650	> 16500	> 4.22	6.36	< 2.14	≥4	FAILS TEST
1%	5 min	10^0	> 165	> 165	> 1650	> 16500	> 4.22	6.36	< 2.14	≥4	FAILS TEST
		10^{-1}	> 165	> 165	> 1650	> 16500	> 4.22	6.36	< 2.14	≥4	FAILS TEST

CONCLUSION

SUSPENSION TEST ACCORDING TO EN 13624:2013

Chemical disinfectants and antiseptics - Quantitative suspension test for the evaluation of yeasticidal activity in the medical area - Test method and requirements
(Phase 2 step 1)

TEST PRODUCT IDENTIFICATION

PRODUCT NAME	:	Alco X Quat
ACTIVE SUBSTANCES	:	Ethyl Alcohol Denat. 59% w/v (73,75% v/v) Benzalkonium Chloride (50%): 0.8%
APPEARANCE OF THE PRODUCT	:	Liquid
STORAGE CONDITIONS	:	Room temperature, darkness
TEST CONDITIONS	:	Test conducted at 20°C ± 1 °C
LOT	:	Not provided
METHOD	:	EN 13624:2013
CONTACT TIME	:	5 minutes
DILUTIONS	:	Undiluted (80%), 50%, 1%
STUDY SPONSOR	:	Σ. KOYZOYNAS & ΣΙΑ Ε.Ε. - ALCOFARM MEDICAL
PRODUCT SUPPLIER	:	Σ. KOYZOYNAS & ΣΙΑ Ε.Ε. - ALCOFARM MEDICAL
PRODUCT MANUFACTURER	:	Σ. KOYZOYNAS & ΣΙΑ Ε.Ε. - ALCOFARM MEDICAL
RECEIPT DATE	:	10/07/2020
STUDY PERIOD	:	11/08/2020-13/08/2020
LAB ID	:	2020-7622/20 23 00657

METHODOLOGY ABSTRACT


A sample of the product as delivered and/or diluted with hard water (or water for ready to use products) is added to a test suspension of fungi in a solution of an interfering substance. The mixture is maintained at 20 °C for 5 minutes. At the end of this contact time, an aliquot is taken, and the fungicidal and/or yeasticidal activity in this portion is immediately neutralized or suppressed. The numbers of surviving fungi in each sample are determined and the reduction is calculated.

RESULT

The product under test: "Alco X Quat" demonstrated yeasticidal activity for surface disinfection (> 4 log reduction) according to EN 13624:2013, at 20±1 °C, under clean conditions for product dilution:

Undiluted (80%) for 5 minutes contact time using as test organism the reference strain: *Candida albicans*.

For the QACS Ltd Laboratory,

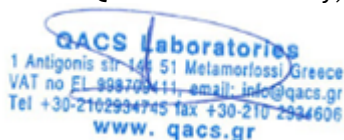

Signature date: 15/09/2020

Lagiopoulos Giorgos

Agronomist-Food Technologist M.Sc.

Study Manager

STUDY SUMMARY / ABSTRACT

SUSPENSION TEST ACCORDING TO EN 13624:2013

Chemical disinfectants and antiseptics - Quantitative suspension test for the evaluation of yeasticidal activity in the medical area - Test method and requirements
(Phase 2 step 1)

TEST PRODUCT IDENTIFICATION

PRODUCT NAME : Alco X Quat
ACTIVE SUBSTANCES : Ethyl Alcohol Denat. 59% w/v (73,75% v/v)
Benzalkonium Chloride (50%): 0.8%
APPEARANCE OF THE PRODUCT : Liquid
STORAGE CONDITIONS : Room temperature, darkness
TEST CONDITIONS : Test conducted at 20°C ± 1 °C
LOT : Not provided
METHOD : EN 13624:2013
CONTACT TIME : 5 minutes
DILUTIONS : Undiluted (80%), 50%, 1%
STUDY SPONSOR : Σ. ΚΟΥΖΟΥΝΑΣ & ΣΙΑ Ε.Ε. - ALCOFARM
MEDICAL
PRODUCT SUPPLIER : Σ. ΚΟΥΖΟΥΝΑΣ & ΣΙΑ Ε.Ε. - ALCOFARM MEDICAL
PRODUCT MANUFACTURER : Σ. ΚΟΥΖΟΥΝΑΣ & ΣΙΑ Ε.Ε. - ALCOFARM MEDICAL
RECEIPT DATE : 10/07/2020
STUDY PERIOD : 11/08/2020-13/08/2020
LAB ID : 2020-7622/20 23 00657

TEST MICROORGANISMS

Candida albicans ATCC 10231
Aspergillus brasiliensis ATCC 16404

RESULT

The product under test: "Alco X Quat" demonstrated yeasticidal activity for surface disinfection (> 4 log reduction) according to EN 13624:2013, at 20±1 °C, under clean conditions for product dilution:

Undiluted (80%) for 5 minutes contact time using as test organism the reference strain: *Candida albicans*.

Results refer to the sample as received and analyzed on the period specified above.

The test report shall not be reproduced except in full, without written approval of the laboratory.

The samples will be stored by the laboratory during 1 month from the end test date.

The study report and raw data will be stored by the laboratory for 5 years.

End of Test Report